


Thank you for purchasing the Sonic-FX Sonic Toothbrush. This toothbrush is equipped with a UV sterilization device. **Please read the entire manual before using and use in accordance with this instruction manual. Keep this manual in a safe place.**




**Warning: Ignoring these warnings and using the product incorrectly may void the manufacturer warranty and/or cause serious injury, illness, or death.**


- Store product out of reach of children.
- Do not put product in or near fire or high temperatures.
- Do not look at the UV lamp. Looking at the UV lamp may have adverse effects.
- Always ensure sanitizer cover is properly closed when in use.
- Do not use product under adverse conditions such as when there is smoke or a burning odor as it may lead to electrical shock or cause a fire.
- If fire or smoke occurs, turn off the toothbrush and unplug the unit from power outlet immediately.
- Do not put a heavy object on the power cord and do not damage or forcibly pull, bend, or twist the power cord.
- If the power cord is damaged, discontinue use immediately or electrical shock or fire may occur.
- Do not alter, damage, or attempt to repair the product. To do so may lead to product malfunction, electrical shock, or injury.
- This product is intended for use by individuals in a general household. Do not use for commercial, professional or any other purposes.
- The charging base is not waterproof. Do not drop or immerse in water. To do so may cause malfunction or electrical shock.
- Do not scratch the UV lamp or rub it with hard or abrasive materials.
- Do not rub the rubber parts of the buttons with hard or abrasive materials.
- Do not store under high temperatures or expose to direct sunlight.
- Do not place on unstable surfaces.
- Clean the toothbrush with water, remove from contact with water, and then place in sterilization device.
- Clean the base of the product body regularly.
- Do not use the same brush for more than one person.
- Do not use in combination with abrasive prophylaxis pastes or toothpastes.
- Do not brush strongly on teeth or gums and do not brush one area for long periods of time or it may cause harm to enamel or gums or deform the brush.
- If you experience bleeding after one week of using this product, discontinue use and consult your dentist.
- Consult a dentist before using this product or stop using the product if any of the following apply:
  - If you have sensitive teeth or gums
  - If you had recent oral surgery or any surgery in or around the mouth
  - If you have bleeding gums or any bleeding in the mouth
  - If you are undergoing dental treatment such as dentures, implants, orthodontics, etc.
- Do not connect or disconnect the power adapter with your bare hands or you may experience electrical shock.
- Do not unplug the power adapter by pulling on the power cord. This may damage the power cord and/or may cause electric shock or fire.
- Use only with the supplied power adapter or malfunction or fire may occur.
- Unplug the power adapter if the product will not be in use for a long period of time or malfunction or fire may occur.
- Only use with the supplied charger or malfunction may occur.

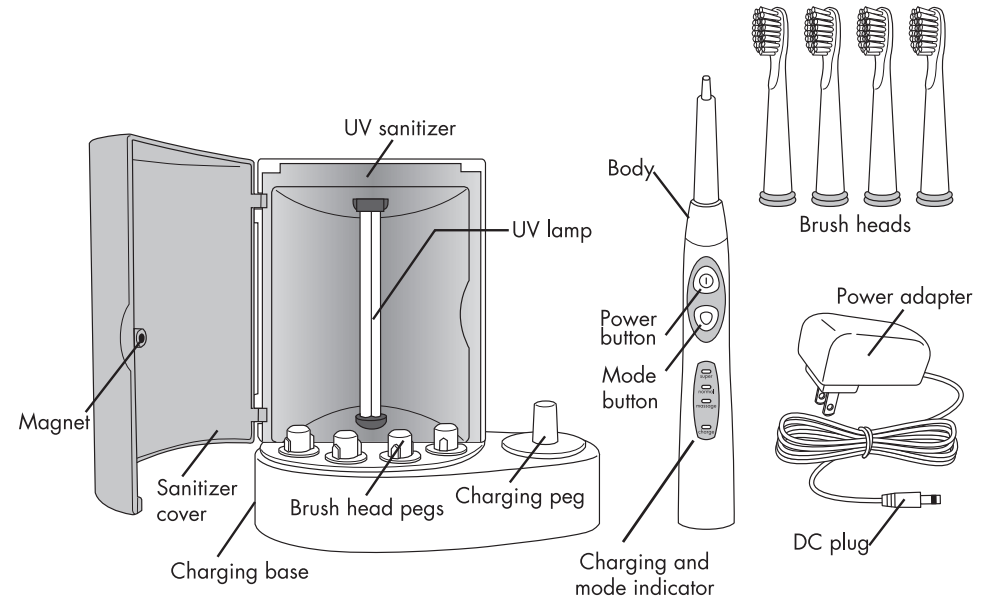

**Install Brush Head:** To install brush head, push brush head all the way down over the tip of the toothbrush body as shown. Give a slight twist to the left until the brush head clicks securely into place. Bristles must be aligned with the power buttons on the brush body. To remove brush head, twist right to release and then remove.

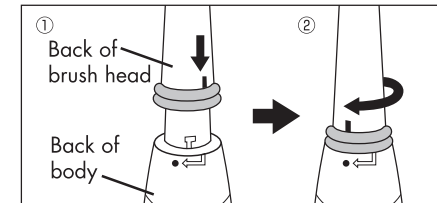

**Charge Toothbrush (Please charge as soon as possible after purchasing):** Connect the power adapter with the back of the charging base as shown and place the unit on a flat surface and in a location where there is no water spray. Plug into power outlet. When the toothbrush is placed on the charging peg, the charging indicator lights up while charging and turns off when fully charged. The base of the toothbrush must be in full contact with charging base in order for charging to take place. **Ensure that toothbrush is charged for 18 hours before first use.** Leave the charge stand plugged into power outlet for normal daily use.

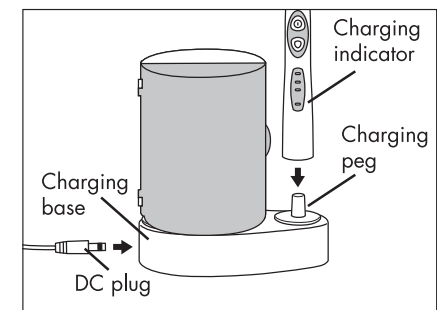


**Warnings:**

- Use only with power adapter supplied with this product or malfunction may occur.
- Base unit is not waterproof. Do NOT immerse in water under any circumstance.

# SONIC-FX™

SONIC TOOTHBRUSH  
WITH UV SANITIZER AND 4 BRUSH HEADS


**To Use Sonic-FX Toothbrush:** Wet brush head with water and add toothpaste. First put the brush head in your mouth and then press power button to start. Press mode button to change modes from Super (high speed) to Normal (medium speed) to Massage (alternating speeds). Start brushing teeth and continue for 2 minutes using the Quadrant Method as shown, in which you brush each quadrant of the mouth for 30 seconds. Sonic-FX will notify you every 30 seconds with a unique pause in vibration after each 30-second interval. Store toothbrush on charging peg between brushings.

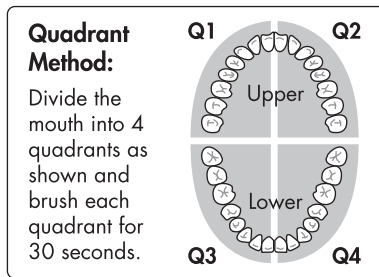

### Warnings:

- It is **not** necessary to exert much pressure or move the toothbrush in the usual tooth brushing motion while brushing teeth with the Sonic-FX toothbrush. Slightly press brush head onto the surface of the teeth and allow the ultrasonic vibrations to clean and polish the teeth while moving the brush head up and down gently over all tooth surfaces, one quadrant at a time. Pressing too hard or brushing in one spot for too long may harm teeth or gums.
- Do not use in combination with tooth abrasives.
- To avoid spray, place brush head in mouth before pressing power button.
- To avoid button malfunction, press buttons with pad of finger only, not fingernails.

**To Sterilize Brush Heads:** Place brush heads on the brush head peg inside the UV sanitizer with the bristles facing the UV lamp. When the cover is closed, the UV lamp will light up and sterilization will start. After 3-5 minutes, the UV lamp will turn off. Sterilization will start every time the case is opened/closed. If lamp does not turn on after closing case, wait 30 seconds and open/close the case again. Lamp should start if everything is properly connected.

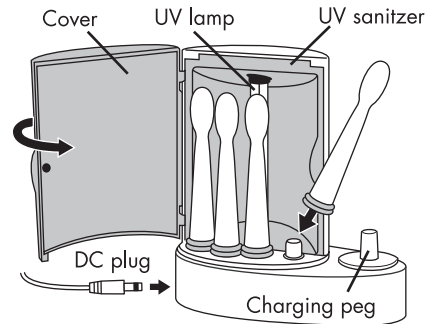

### Warnings:

- Store product out of reach of children.
- Avoid looking at the UV sanitizer while lamp is on.
- Avoid exposure to the UV lamp while lamp is on.
- Do not allow UV lamp to come in contact with water.
- Cover must be properly closed or it will not function correctly.



### Note:

- This product is intended for normal household use. Do not use for commercial, professional or any other use.

### Troubleshooting:

- **Brush does not power on:** Check all connections. If still not working, charge for at least 24 hours.
- **Brush or sanitizer functions strangely:** Battery may be low. Charge for at least 24 hours.
- **Brush does not stay on long enough to brush teeth:** Be sure to store toothbrush on charging peg between brushings. Make sure unit is kept plugged in.
- **UV lamp does not come on:** Make sure UV sanitizer case is closed properly. Door latch is magnetic and should securely snap shut. Check to make sure the door magnet is in place.

### To Clean:

- Rinse bristles after every use. Replace brush heads after approximately 3 months of daily use, or the equivalent thereof, or as needed.
- Clean toothbrush body, UV sanitizer exterior, and charging base with a soft damp cloth on a regular basis to ensure proper working condition.
- Do not use volatile cleaning solvents.
- Do not submerge charging base/UV sanitizer in water.

### Specifications:

Materials	ABS (body), nylon (brush)	Time Needed for Sterilization	3-5 minutes; automatic shut-off
Body Size	About 126Wx62Dx258H mm (w/charging seat)	Number of Brushes in the Sterilizer	Up to 4 brush heads
Body Weight	About 250g (excluding the adapter)	Speed Adjustment	3 stages: Super, Normal, and Massage
Rechargeable Battery	Lithium-ion	Hours of Continuous Use	About 14 days on full charge; approx. 2 min. of use 2x/day
Charging Mode	Inductive charging	Power Consumption	About 5W when charging; about 35W when sterilizing
Charging Time	Approximately 10-18 hours	Accessories	3 replacement brush heads

### Warranty

Product Name: Sonic-FX Sonic Toothbrush

Warranty Period: Your Sonic-FX Toothbrush comes with a 90-day warranty that begins on the date of purchase.

### Terms:

- 1) During the warranty period, if the product malfunctions during normal use in accordance with this instruction manual, the product will be repaired or replaced free of charge.
- 2) The product will not be repaired or replaced free of charge even during the warranty period if any of the following occurred:
  - a) malfunction and/or damage due to improper use, repair, or reconstruction
  - b) malfunction and/or damage due to improper operations, fall, or improper shipment after purchase
  - c) malfunction and/or damage due to natural disaster such as fire, earthquake, flood, etc., or abnormal voltage
  - d) malfunction and/or damage due to the malfunction and/or damage of other devices connected to this product
  - e) malfunction and/or damage due to use for other than normal personal use.
- 3) This warranty is valid only in the USA by contacting Teeth Whitening Technologies at [www.TeethWhiteningTechnologies.com](http://www.TeethWhiteningTechnologies.com). We must be able to locate your purchase in our database or you may be required to provide proof of purchase including date of purchase.

### Technical Assistance:

1205 NE 95th St. STE A  
Vancouver, WA 98665

To contact us or order replacement brush heads, please visit:

[www.SonicFXToothbrush.com](http://www.SonicFXToothbrush.com)